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
This device is a smart pedal exerciser, which can be folded and easily stored. Its use is recommended to promote blood circulation and support muscle development. It is a fun and enjoyable way to perform cardiovascular exercises. Placed on a table or directly on the floor, it can be used with the arms or legs..

DESCRIPTION

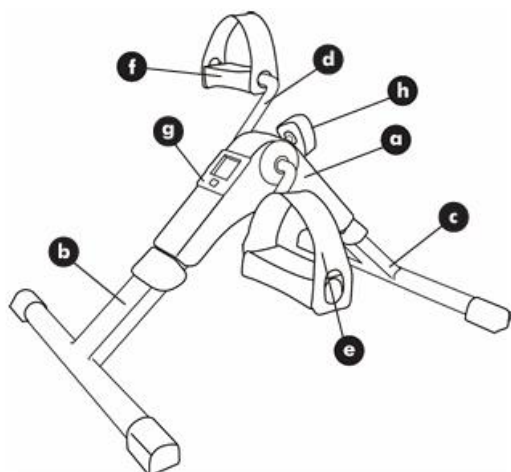- a- Main body
- b- Front base
- c- Rear base
- d- Handlebar
- e- Right pedal
- f- Left pedal
- g- Counter and battery compartment
- h- Resistance regulator

ASSEMBLY

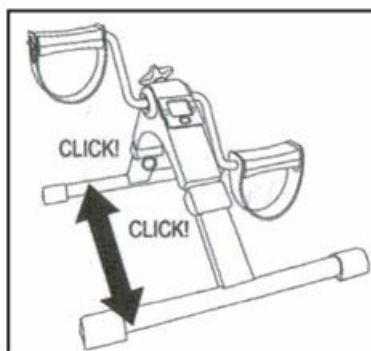
fig. 1

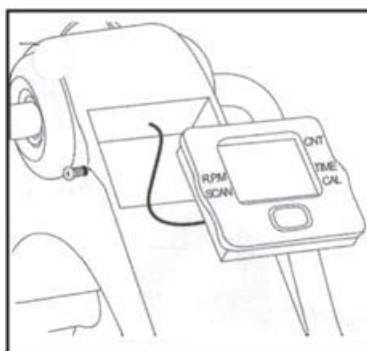
fig. 2

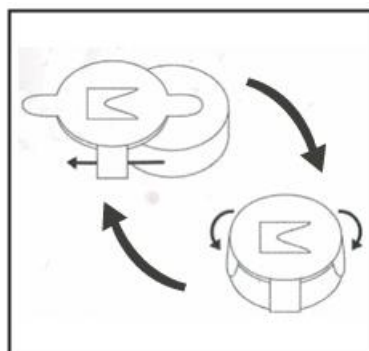
fig. 3

1. Unfold the front and rear bases until you hear a "click." This indicates that the structure is properly secured for exercising. Before starting any exercise, check that both bases are well secured (see fig. 1).
2. Insert the display unit into its slot. Make sure the wire is correctly positioned inside (see fig. 2).
3. Place the intensity regulator screw in its designated slot, located on the upper part of the mini exercise bike.

SAFETY PRECAUTIONS

- These instructions are very important. Please read this instruction manual carefully before using the device and keep it handy for future reference.
- IDENTITES is not responsible for any misuse of the product or any use other than what is described in this manual.
- Before using the product for the first time, unpack it and check that it is in perfect condition. If not, do not use the device and contact IDENTITES, as the warranty covers any manufacturing defects or original damage.
- This product is intended for home use only. Do not use the device outdoors. Never expose the device to weather conditions such as rain, snow, or direct sunlight.
- This device may be used by children aged 8 and older, as well as disabled individuals or persons with reduced sensory or mental capabilities, or those with limited experience and knowledge, provided they are properly supervised and instructed and understand the risks involved.
- Cleaning and maintenance operations that can be performed by the user should not be carried out by children, unless they are over 8 years old and under supervision.
- Keep the device out of reach of children under 8 years old.
- Children must be supervised to ensure they do not play with the device.
- Do not immerse the device in water or any other liquid. Do not mix new batteries with used or discharged batteries.

USE

This pedal exerciser has four different functions:

- **CNT:** Total number of rotations performed during this session.
- **TIME:** Effective exercise time during this session.
- **CAL:** Calories burned during the workout.
- **RPM:** Displays the average rotations per minute based on exercise intensity.

Before starting the exercise:

Place the device on a flat surface and adjust the pedaling resistance according to the desired intensity (h). To increase resistance, turn the knob counterclockwise.

Usage:

1. Powering on:

Press the red button located on the display (g) or start pedaling. The display (g) will turn on automatically.

- The screen displays "STOP" when you stop pedaling and disappears as soon as you resume exercise.
- The display automatically turns off a few minutes after stopping pedaling.

2. Viewing results:

- **Automatic mode:**

Press the red button once. The word "SCAN" blinks, and all the results (TIME, RPM, CALORIES) are displayed sequentially, with each value appearing for about 5 seconds.

- **Manual mode:**

To select a specific function to display during the exercise, press the red button twice to choose the desired function (TIME, CAL, RPM). If you press the red button once, the device returns to automatic mode without selecting a specific function.

- **Resetting the display:**

To reset the values to zero, hold the red button for a few seconds. All indicators will reset to "0".

3. Replacing the batteries:

Remove the display unit (g) from its original position using the slots around the display. This will give access to the battery compartment, which can be replaced without using tools (see fig. 3).

CLEANING AND STORAGE

- First, clean with a damp cloth, then dry with a soft cloth.
- Store the device in a safe place.
- To fold and store the device, pull the rings until both bases detach.

TECHNICAL SPECIFICATIONS

- **Required battery:** 1 x LR44, 1.5V
- **Maximum load capacity:** 90 kg

This device complies with safety requirements and the provisions of Directives 2014/35/EU on low-voltage appliances and 2014/30/EU on electromagnetic compatibility.

RECYCLING

Packaging

The materials used for this product's packaging are part of a waste collection, sorting, and recycling system. If you wish to dispose of it, please use the appropriate bins or municipal collection services designated for each type of material.

Electrical and Electronic product

This symbol indicates that if you wish to dispose of this product at the end of its life, you must take it to an authorized collection facility for the separate disposal of Waste Electrical and Electronic Equipment (WEEE).

This product does not contain concentrations of materials considered harmful to the environment.

WARRANTY

This product is guaranteed for 2 years from the date of purchase by the user. IDENTITIES will provide replacement parts for any part that has broken due to a manufacturing, design or material defect. Contact the distributor from whom you purchased the product in case you need to make a claim. Failures resulting from accidental damage, misuse, alteration, or wear and tear are not covered by this warranty. Contact the dealer from whom you purchased the product to apply the warranty.

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