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INSTRUCTIONS MANUAL EN

IMPORTANT

Read the instructions carefully before using this product and keep them if necessary.

User: If you are unable to read or understand any warnings, precautions or instructions, please contact your health care professional or distributor before using this product or risk and damage may result.


Consult your doctor or health care professional before using a walking aid. He/she will advise and inform you on the most appropriate walking technique for you.

Distributor: this leaflet must be given to the users of this product.

PRODUCT DESTINATION

Tasmania bed tables 2 tops are medical devices. They allow the user to read or take a meal in bed or in an armchair. They are adjustable in height and tilt to adapt to the needs of the user during the day.

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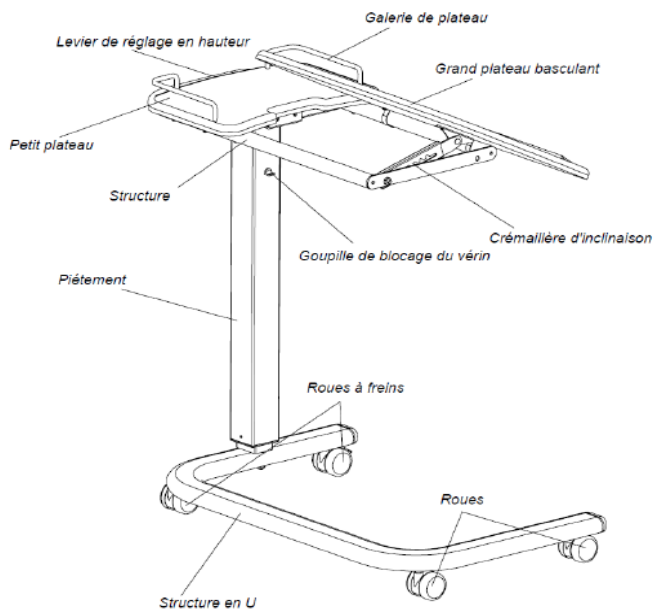SAFETY PRECAUTIONS


Failure to follow the safety precautions may result in the risk of serious injury.

- Use the bedside table only as described in this manual.
- Install the bed table on a flat surface.
- Store the bedside table in an indoor environment, away from moisture and dust.
- This bedside table uses a high-pressure pneumatic cylinder to adjust the height of the bed top. Do not remove the locking pin from the column assembly before use.
- Before use, make sure that all parts are present and in good condition. In case of damage, please contact your dealer.
- Carefully follow the assembly instructions and ensure that all components are properly secured before use.
- Before use, tighten the brakes.
- In case of damaged, broken or worn parts, please contact your dealer

MOUNTING


INSTALLATION

Place the bedside table on a clear surface and remove the packaging. Do not discard the packaging before the assembly is completed.

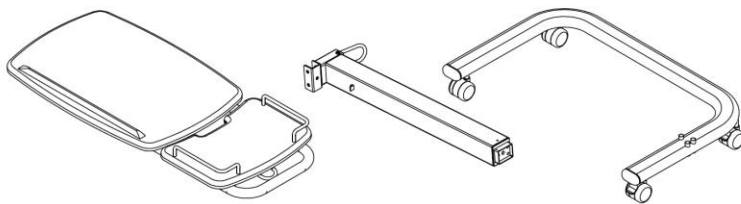
STEP 1

Place the table upside down on a flat surface as shown. Position the frame above the table top as indicated. Insert the hexagon socket screws through the frame into the table top frame as indicated. Use the M8 hexagonal key to tighten the screws.

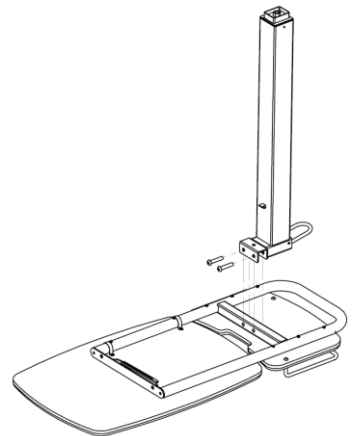STEP 2

Place the U-shaped frame upside down over the base as indicated on the left. Screw the 3 screws through the plate, the base and the end of the base.

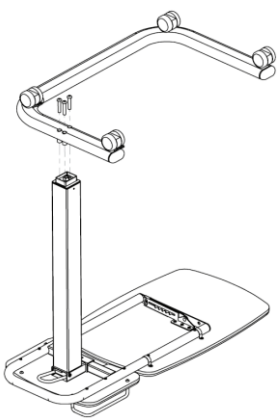STEP 3

Turn the table over and place it on the floor. Make sure that all parts are perfectly tightened before using the bed table.

Remove the locking pin from the frame.




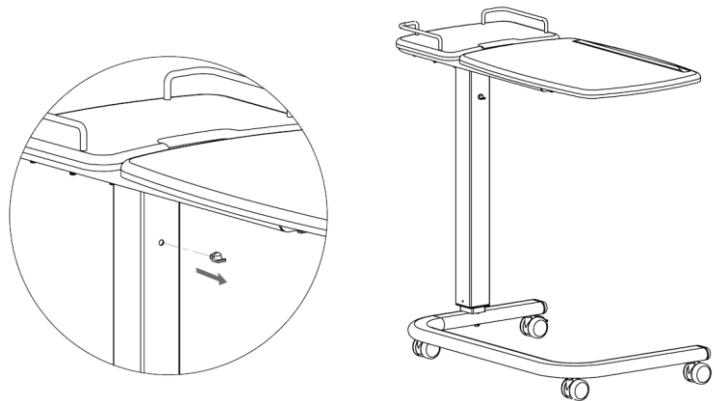
TABLE SLOPE

The Tasmania 2-top bedside table must be placed on the floor before tilting the tabletop.

1. Pull the large tray out to tilt it. The ratchet mechanism allows 6 tilting positions.
2. Tilt the large tray to the desired angle.
3. To return the large tray to a flat position, tilt it to its maximum position to release the ratchet and then gently place the large tray back in place.

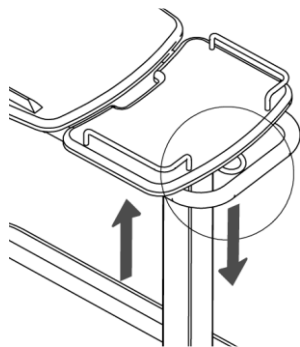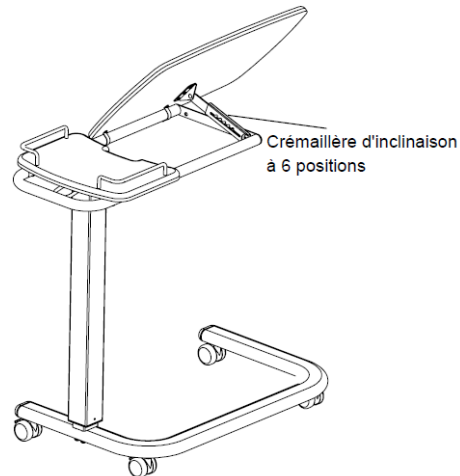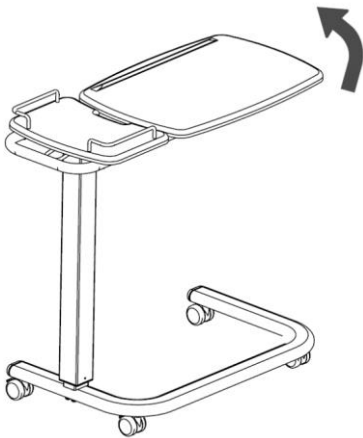HEIGHT ADJUSTMENT

Simply pull the table top upwards to adjust the table height.

To lower the table, pull the adjustment lever, adjust to the correct height and then release the lever.

CARE AND MAINTENANCE

Failure to care for and monitor your product may result in the risk of serious injury.

- Clean the product regularly with lukewarm water, a soft cloth and a neutral detergent.
- Do not use cleaning products containing phenol or chlorine, as this may damage aluminium and polyamide materials.
- Dry with a clean, dry soft cloth.

- Make sure there is no accumulation of dirt or debris around the axle, brakes or wheels.
- Periodically check all screws and bolts for wear and strength. Replace all broken, damaged or worn parts immediately.
- Periodically check the cylinder.
- In the event of a malfunction, do not repair the table yourself. Consult an authorised dealer.

CHANGE OF OWNERSHIP

The product can be reused. It must first be cleaned and disinfected in accordance with hygiene instructions. When transferring the product, remember to hand over all the necessary technical documents to the new user. Beforehand, the product must be inspected by an authorised specialist.

RECYCLING

When the product has become unusable and you have to throw it away: please contact your specialist dealer. If you wish to carry out recycling for yourself, ask your local waste management company for disposal guidelines.

STORAGE

Failure to observe the storage conditions may lead to deterioration of the product and thus to the risk of serious injury.

Do not store your product for an extended period of time near a heat source or in the sun (e.g. behind a window or near a radiator) or near a cold source.

Keep away from all flames and sources of sparks.

Observe the conditions for storing and stowing your product.

- In a dry and temperate place
- Protect your product by packaging it from dust, corrosion (e.g. abrasive elements, sand, sea water, salty air).
- Store all removed parts together in the same place (or locate them if necessary) to avoid mixing them with other products when reassembling.
- All components must be stored without any load (do not place any objects that are too heavy on the product elements, do not jam anything...).

TECHNICAL SPECIFICATIONS

Base dimensions	Width 71 x Depth 45 cm
Large tray dimensions	57 cm x 39 cm
Dimensions of small tray	25 cm x 35 cm
Adjustable height	79 cm à 114 cm
Clearance between U-shaped frame	59 cm
Total weight	16 kg
Maximum supported weight	25 kg

PERIOD OF USE

The life expectancy of this product is 5 years under normal use, safety and maintenance conditions. Beyond this period, the product can be used for as long as it is still in good condition.

WARRANTY

This product is guaranteed for two years from the date of purchase by the user. IDENTITIES will provide replacement parts for any part that has broken due to a manufacturing, design or material defect. Contact the distributor from whom you purchased the product in case you need to make a claim. Failures resulting from accidental damage, misuse, alteration, or wear and tear are not covered by this warranty. Contact the dealer from whom you purchased the product to apply the warranty.

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