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EN

IMPORTANT

Read the instructions carefully before using this product and keep them if necessary.

User: If you are unable to read or understand any warnings, precautions or instructions, please contact your health care professional or distributor before using this product or risk and damage may result.

Consult your physician or health care professional before using any assistive device. He/she will advise and inform you on the most appropriate use in your case.

Distributor: this user manual must be given to the users of this product.

PRODUCT DESTINATION

The Magnus rollator is a walking aid designed for people with reduced mobility. This rollator is intended for outdoor use and can be folded for easy storage or transport.

SAFETY PRECAUTIONS


Failure to follow the safety precautions may result in the risk of serious injury.

- Do not exceed the maximum weight limit of 150 kg.
- Do not exceed the maximum load limit of the basket (5 kg). Loads must only be carried in the basket. Do not hang heavy objects on the rollator, as this may affect its stability.
- Ensure that all parts (wheels, cables, handles) are securely fastened and in good working condition.
- Do not use the rollator on stairs or escalators.
- Use only accessories and parts supplied by the manufacturer.
- The rollator is a walking aid. **It is not a wheelchair.** It is not designed to transport a seated person. **Such use may damage the wheels and the structure of the rollator**, weakening it and preventing it from meeting its mechanical specifications in future use. Use the rollator only as a walking support.
- Do not use the rollator to walk backwards.
- Do not make any adjustments to the rollator while it is in use.
- Check that the brakes are functioning properly before use.
- All wheels must always be in contact with the ground.
- When the rollator is **stationary**, the brakes must be locked. Do not attempt to move the rollator if the brakes are locked.

- Make sure to place your hands on each handle so that the rollator remains in a stable, parallel position.
- Take care when going up curbs or when using inclined, uneven, or slippery surfaces. Failure to follow this instruction may result in a significant risk of falls or injury.
- Do not use the rollator on wet, icy, or snowy slopes.
- In case of damaged, broken, or worn parts, contact your dealer.
- Check that the folding mechanism is properly secured and that the height adjustment handles are locked before use.
- The device must not be used by persons with severe balance disorders or whose muscle strength is insufficient to stand and walk, even with the assistance of a walker. The need for and ability to use a walker must always be assessed and validated by a doctor or physiotherapist.
- Ensure that this user manual is read by all people using the device. The manufacturer declines all responsibility for damage and/or injury resulting from failure to follow the instructions in the manual.
- The brakes **MUST** be locked **BEFORE** using the seat.
- Use the device only for its intended purpose.
- Do not make any structural modifications to the device without prior written consent from the manufacturer.
- When adjusting the device, the user must take care not to pinch fingers or other body parts between moving components.
- The device must only be used on a stable surface.
- Beware of the risk of pinching when folding and unfolding the device.

PRODUCT DESCRIPTION

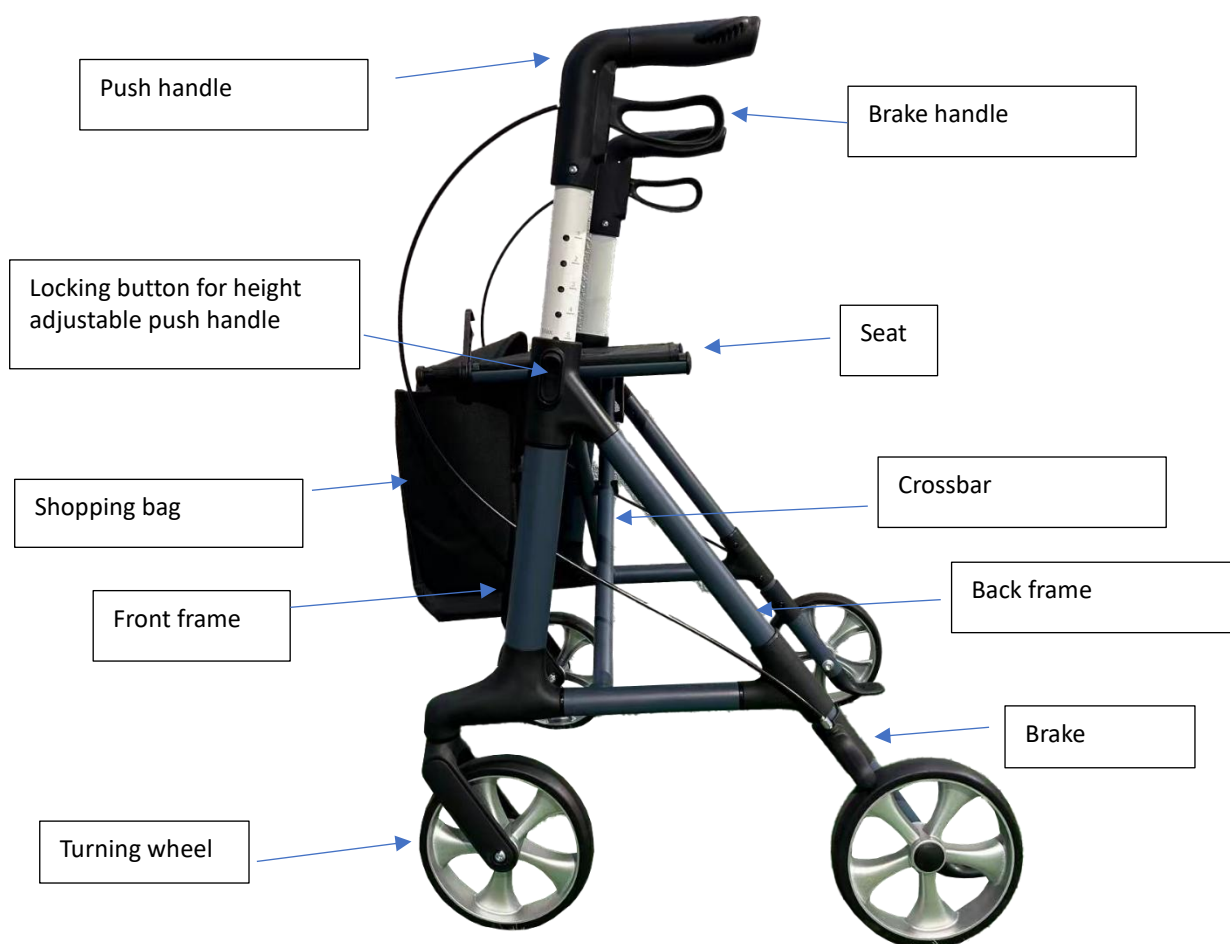ASSEMBLY & USE

With reference to the illustration above, the device can be prepared for use intuitively. However, we recommend following the steps below :


Initial check :

- First, remove the device from its packaging and check (using the illustration above) that all parts are present.
- If you notice that a part is missing or damaged, do not assemble the product. Contact your dealer immediately.


Complete package contents: One seat, one basket, one user manual, and one main support mechanism.

ASSEMBLY :

Step 1 :

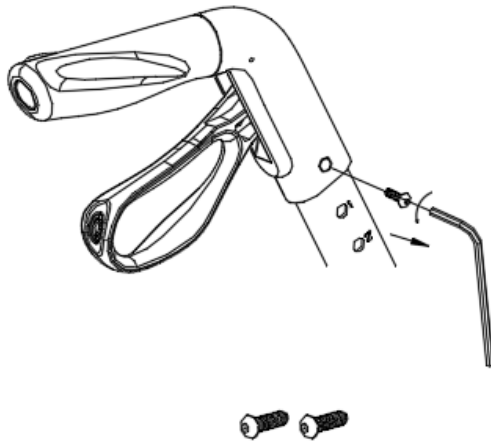Step 2 :

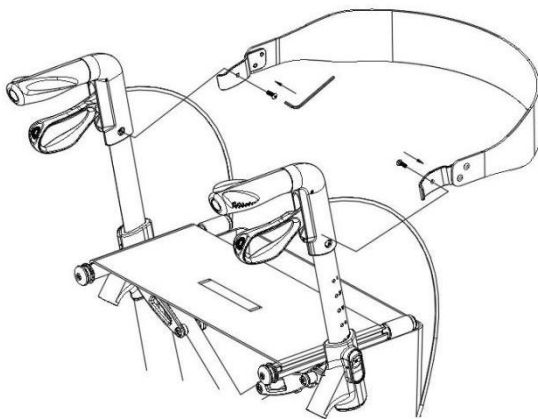Opening the rollator:

1. Open the locking latch of the folding system.
2. Hold the handles with both hands and pull them apart on each side.
3. Press down on the support tubes on each side until you hear a “click,” indicating that the mechanism is properly locked.

- After unfolding or assembling the rollator, always ensure that the seat support is properly locked and that the frame no longer moves.
- Make sure that both push tubes are securely locked into the main support tube using the height adjustment locking button, to prevent any accidental slipping during use.
- Always check that the accessories are properly and securely attached BEFORE use.
- Failure to follow these instructions may result in injury.

Folding the rollator :

1. Adjust the support tubes on both sides to the lowest position.
2. Pull the folding latch located in the center of the seat until the left and right sides of the frame are fully folded.
3. Secure the seat support tubes on both sides using the folding system locking latch.

Handle height adjustment :


- Unlock the handle adjustment knob.
- Adjust the height as needed.

- Both handles must be set at the same height.
- The correct handle height should correspond to the user's hip level.
- Never leave the device without activating the brakes.

Using the hand brakes :

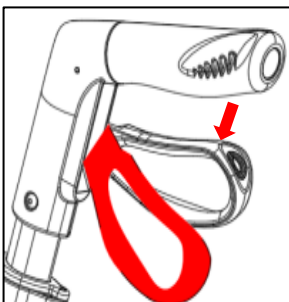
1. To stop the rollator: simply pull the brake handles upward, toward the handles.
2. To remain stationary: push the brake handles down to lock them in place.
3. To move again: release the brake handles. They will automatically return to their original position and the brakes will disengage.

- Always activate the brakes when you do not intend to move the device.
- Never leave the device without the brakes engaged when someone is sitting on it.
- Do not use the rollator on a slope without engaging the brakes.
- When seated, sit upright with your forearms or hands resting lightly on the handles.
- Never lean backward when seated on the rollator.

Locking / Unlocking the hand brakes:

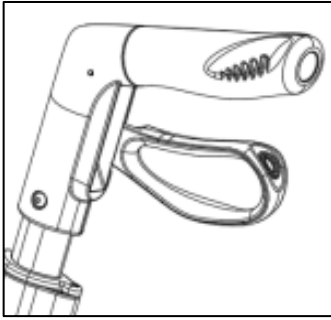1. To lock the brakes, push the brake handles downward until you hear a “click.”
2. To unlock the brakes, pull the brake handles upward until they release.

When using the rollator in a stationary position, the hand brakes must always be locked.

CARE AND MAINTENANCE


Failure to care for and monitor your product may result in the risk of serious injury.

- Regularly clean the rollator with lukewarm water and a soft cloth. Do not use abrasive products, detergents, or solvent-based cleaners.
- Dry with a clean, soft, dry cloth.
- Ensure that there is no accumulation of dirt or debris around the axle, brakes, or wheels.
- Do not use cleaning products containing phenol or chlorine, as this may damage the rollator.
- Lubricate the wheels with a small amount of oil if they squeak.
- Periodically check the wear and tightness of all screws and bolts. Replace any broken, damaged, or worn parts immediately.
- Check that the brakes are functioning before each use.
- Check the parking brakes by pushing the handles into the locked position, then attempting to push the rollator.
- Ensure that all parts of the rollator are securely fastened at all times.
- In case of malfunction, do not attempt to repair the rollator yourself. Consult an authorized dealer.

Periodic safety inspection:

To ensure safe use, please check the following points for any loose bolts or nuts before and after each use.

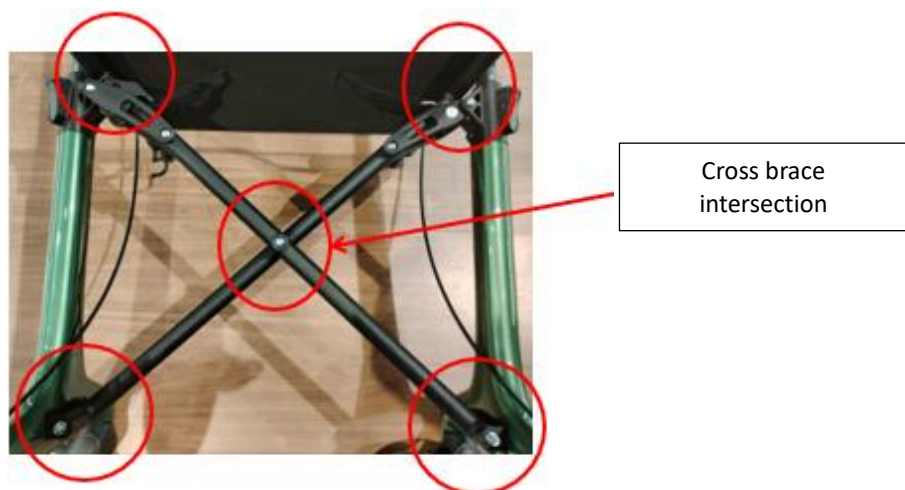1. Cross brace intersection

Check the bolt, washers (including the saddle washer), and nut at the center of the X-shaped frame to ensure they are not loose.

2. Upper connections

Check the bolts, flat washers, and nuts at the connection points between the upper part of the cross brace and the rollator frame on both the left and right sides for any looseness.

3. Lower connections

Check the bolts, flat washers, and nuts at the connection points between the lower part of the cross brace and the rollator frame on both the left and right sides for any looseness.

Safety warning

If any fasteners are loose, tighten them immediately using appropriate tools (not supplied) before using the rollator, in order to avoid any potential risk.

- **Checklist :**

To keep the rollator in good condition, check before each use, and check regularly as follows:

Item to check	At all times	Weekly	Monthly
Brakes (operation, locking, cables)	X		
Wheels (rotation, wear, fastening)		X	
Seat (stability, condition, fastening)		X	
Backrest (fastening, support)		X	
Height adjustment push buttons (proper locking, left/right symmetry)		X	
Folding system (full frame locking)	X		
Structural welds (especially around the wheels)			X
Handles (fastening, stability, safety)	X		
Screws and fastening elements (general)		X	
Overall frame stability			X

CHANGE OF OWNERSHIP

The product can be reused. It must first be cleaned and disinfected in accordance with hygiene instructions. When transferring your product, remember to hand over all the necessary technical documents to the new user. Beforehand, the product must be inspected by an authorized specialist.

RECYCLING

When the product has become unusable and you have to throw it away: please contact your specialist dealer. If you wish to carry out recycling for yourself, ask your local waste management company for disposal guidelines.

STORAGE

Failure to observe the storage conditions can lead to deterioration of the product and thus to the risk of serious injury.

- Do not store your product for an extended period of time near a heat source or in the sun (e.g. behind a window or near a radiator) or near a cold source.
- Keep away from all flames and sources of sparks.
- Respect the conditions for storing your product.
- In a dry and temperate place.
- Protect your product from dust, corrosion (e.g. abrasive elements, sand, sea water, salty air).
- Store all removed parts together in the same place (or locate them if necessary) to avoid mixing them with other products when reassembling.
- All components must be stored without any load (do not place any objects that are too heavy on the product elements, do not jam anything...).

PERIOD OF USE

Life expectancy of this product is 5 years under normal use, safety and maintenance conditions for a user. Beyond this period, the product can be used for as long as it is still in good condition.

If the device is used in a multiuser environment (such as a nursing home or hospital) then service life may be reduced.

TECHNICAL FEATURES

Overall dimensions	W. 62 × D. 65 × H. 78 to 94 cm
Seat dimensions	W. 44.5 × D. 23 × H. 60 cm
Maximum bag load	5 kg
Weight	6,8 kg
Maximum user weight	150 kg
Folded overall height (minimum)	88,8 cm
Folded width	22 cm
Wheels (diameter × width)	Black PA + 15% GF + TPR SOFT surface

GUARANTEE

This product is guaranteed for two years from the date of purchase by the user. IDENTITIES will provide replacement parts for any component that fails due to a defect in materials, manufacturing, or design. Contact the distributor from whom you purchased the product in case you need to make a claim. Failures resulting from accidental damage, misuse, alteration, or wear and tear are not covered by this warranty. Contact the dealer from whom you purchased the product to apply for the warranty.

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